English	5
Bahasa Melayu	7
Čeština	10
Crnogorski	12
Dansk	14
Deutsch	16
Eesti	19
Español	21
Español	24
Français	26
Français (Canada)	29
Hrvatski	32
Italiano	34
Latviešu	37
Lietuviškai	39
Magyar	41
Nederlands	44
Norsk	47
Polski	49
Português	52
Română	54
Shqip	57
Slovenščina	59
Slovensky	61
Srpski	63
Suomi	65
Svenska	67
Türkçe	70

## Battery removal instructions for professionals - Sonicare Prestige 9900

The built-in rechargeable battery must only be removed by a qualified professional when the appliance is discarded. Before removing the battery, make sure that the battery is completely empty.

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.**


**When you handle batteries, make sure that your hands, the product and the batteries are dry.**

**To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.**

To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare toothbrush.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.

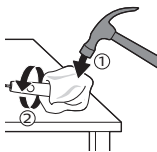
- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

Note: You may have to hit on the end several times to break the internal snap connections.

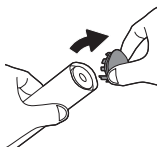
- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

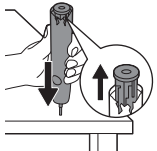
- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

- 6 Pull the internal assembly to detach it from the housing tube.

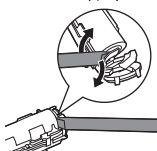- 7 Insert the screwdriver between the bottom of the battery and the frame to break the metal tab connecting the battery to the printed circuit board.

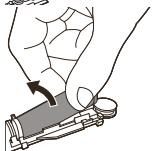- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge.

- 10 The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

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©2021 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.
Sonicare et le logo Sonicare sont des marques de commerce de Philips Oral Healthcare, LLC et (ou) KPNV.
Philips Oral Healthcare, LLC.