



PHILIPS

InnoSpire Go

Portable
mesh nebulizer

Reduces treatment time by 25%*

The easy-to-use and virtually silent InnoSpire Go has been specifically designed to shorten nebulizer treatment times, giving your patients more time for life.

Utilizing established, clinically proven vibrating mesh technology adopted by hospitals worldwide, InnoSpire Go delivers all commonly prescribed respiratory medications.



As easy as 1, 2, 3...



1. Fill

Lift the green lid, pour in the medication and the InnoSpire Go is ready to use.



2. Treat

Press the on button and the InnoSpire Go will begin to nebulize.



3. Clean

Remove the mouthpiece and quickly rinse it in warm water.

* Versus predecessor Aeroneb Go.

InnoSpire Go provides an easy therapy experience, which helps patients administer therapy correctly.



Treatment in 4 mins**

Efficient aerosol therapy with audible and visible indicators that signal therapy completion - before automatically turning off.



Virtually silent

A small, portable and discreet therapy device, allowing your patients to take their therapy when and where they need it.



Built-in battery

Rechargeable, long life, built-in battery that provides up to 30 treatments between each charge.



Simple, only two parts

The easily detachable mouthpiece makes the InnoSpire Go simple to use, clean and maintain, encouraging greater adherence.***

InnoSpire Go aerosol data

	Salbutamol 2.5 ml	Ipratropium bromide 2 ml	Budesonide 2 ml
MMAD μm	3.99	3.93	4.84
Respirable fraction % particles < 5 μm	64.4%	65.3%	52.1%
Treatment time (min:sec)	4.18	3.44	3.34
Residual volume ml	0.3	0.2	0.3




Product specification table

Battery	Lithium ion
Dimensions	7.0 cm x 4.5 cm x 13.5 cm
Weight	111 g
Warranty	Handset 2 years Mouthpiece 1 year

** Using 2.5ml salbutamol

*** A mask is included for children that find it difficult to use the mouthpiece.

Find out more – www.philips.com/respironics

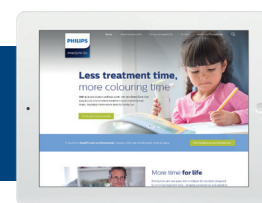

© 2016 Koninklijke Philips N.V. All rights reserved.
Specifications are subject to change without notice.
Respironics Inc. and Respironics Respiratory Drug
Delivery (UK) Ltd are Philips companies.



Respironics Respiratory Drug Delivery (UK) Ltd
Chichester Business Park, City Fields Way
Tangmere, Chichester PO20 2FT UK
KW 09/01/2016 PN452302
PN 1132091 MCI 4107773
www.philips.com/respironics