

PHILIPS

Efficia

CMS200 central
monitoring system



Cost-effective central monitoring

You strive to provide high quality care, on a limited budget, for large numbers of patients. A cost-effective way to centralize your monitoring can make a significant contribution to helping your staff work efficiently. The Efficia CMS200 central monitoring system works with your Efficia CM series or SureSigns VM¹ patient monitors to support not just centralized monitoring and alarming, but also enhanced access to historical patient data. We have designed the Efficia CMS200 with core functionality that allows you to deliver excellent care while tackling cost goals. When budget and quality matter, the choice is Philips Efficia.


Key advantages

- Quick-to-learn, intuitive user interface for ease of use
- Reliability at a cost-effective price, for a low total cost of ownership
- Centralized monitoring, secondary alarming, bi-directional communication of selected settings, and data review including full disclosure data, for efficient workflows
- HL7 output to third party EMR or HIS systems

Flexible and easy to use

The Efficia CMS200 lets you display waveforms, parameters, and alarms from networked patient monitors for up to 32 patients. You can configure the layout of the pages and sectors to make it easy to access the wealth of information at your disposal, whether you are examining the details from up close, or glancing from a distance.

In addition to viewing a patient's physiological data in real-time, you can also review it – including all the waveforms in full disclosure – from up to the previous ten days, for example, to investigate changes in his or her condition.

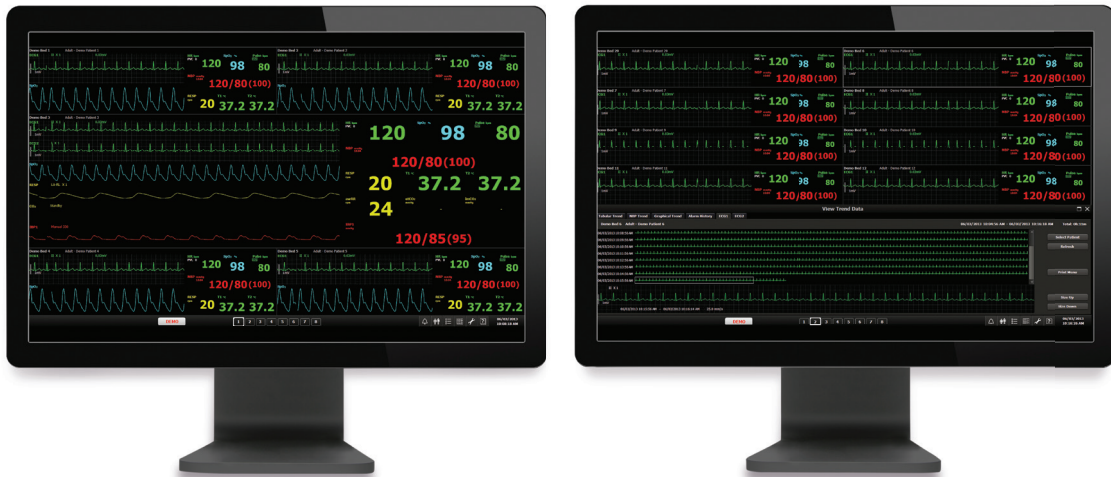

Easy setup and install

The Efficia CMS200 includes everything you need (including keyboard, mouse, speaker, and network switch) to attach to your patient monitoring network and get started right away. The easy-to-understand operation and the integrated Instructions for Use can save you time on training.

Enhance your alarm management capabilities

The Efficia CMS200 is a secondary alarm device, gathering physiological and technical alarms across the network, providing audible and visual notification, and giving you the option of acknowledging them at the central station. In addition, the Efficia CMS200 can forward monitoring data from the networked monitors to your EMR.

Control the monitor with remote access

With the Efficia CMS200 you initiate certain functions on the monitor from the central station. You can share patient demographics, start or stop NBP measurements, and silence alarms from the central. This allows you to focus on providing productive care for your patients from beyond the bedside.

Windows is a registered trademark of Microsoft Corporation in the United States and other countries.

¹ The Efficia CMS200 works with SureSigns VM patient monitors with software revision A.03.93 and later.

