



Philips
Body analysis scale

Body weight measurement
Body fat and BMI analysis
Syncs with Philips health app



DL8780/37

Connected weight tracking

Syncs with free Philips HealthSuite health app

The connected body analysis scale provides accurate weight measurements and shows your body fat percentage and BMI in the free Philips HealthSuite health app so you can conveniently track your body composition over time.

Delivers accurate weight measurements

- Accurate, fast measurement based on 4 high-precision sensors
- BMI and body fat percentage shown in the Philips health app

Personalized tracking in the Philips health app

- Automatically syncs to the Philips health app via Bluetooth®
- Allows up to 8 different users to track their readings
- Monitor body composition over time

Designed to fit your bathroom

- Designed with large screen to clearly show results
- Gives results in kilograms, stones and pounds
- Designed with soft edges for safe use

PHILIPS

Body analysis scale

Body weight measurement Body fat and BMI analysis, Syncs with Philips health app


DL8780/37

Highlights

Fast and accurate


The scale has a simple sense-on technology that starts when you step on it. Four high-precision sensors deliver the accurate weight measurements.

Personalized tracking

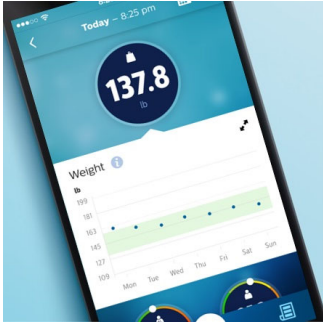

Determines your normal weight over time, along with your body fat percentage and BMI so you can track and spot changes quickly and easily in the Philips HealthSuite health app.

Body composition analysis

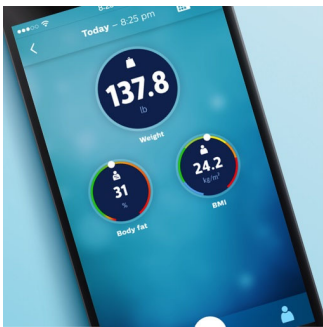

Four high precision sensors deliver accurate weight measurements. Bio-Impedance Analysis (BIA) sends a small unnoticeable current through the body and establishes the body fat percentage. BMI is calculated based on weight and height of the user.

Phone compatibility

Automatically transfer your measurements to the Philips HealthSuite health app on your smart phone, so you can monitor changes over time. Compatible with

iOS8 or higher/Android 4.4 or higher. For an overview of (in)compatible devices please have a look at www.philips.com/personalhealthcompatibility

Multiple user identification

Automatically recognizes you when you step on the scale and sends your measurement data to the Philips HealthSuite health app on your smartphone.

Large screen



Making it easy for you to read your measurement results with a large backlit LCD screen.

Kilograms, stones and pounds

No matter whether you measure in kilograms, stones or pounds, the Philips body analysis scale will give you the results you can understand.

Soft edges



Compact and streamlines, the soft edges provide an appealing look and a safe user experience

Specifications

Syncing

- Phone compatibility: iOS8+(iPhone®4S+)/Android4.4.2
- Data transfer technology: Bluetooth® Smart

Design

- Color: Black
- Dimensions scale: 12.6x12.6x0.98in
- Net. Weight: Approximately 4.37lbs
- Display: LCD display, white backlight
- Display size: 2.91x2.09in
- Soft edges

Ease of use

- Sense-on technology
- Easy switch to kg., st. & lb.
- Large LCD display
- Automatic user recognition
- Can be used as stand-alone

Technical

- Accuracy: $\pm 0.7\text{lb}-110\text{lb}/+1.0\text{lb}$ per 110lb+
- Measurement range: 11 lb to 397 lb
- Measurement units: Kilogram / pound / stone
- Graduation: 0.1lb

Power

- Power source: 6V (4xAAA alkaline batteries)

Accessories

- 4x AAA batteries

Service

- 2-year guarantee



Issue date 2017-04-08

Version: 6.4.1

© 2017 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com